

HORIBA

Scientific

*Industry's Widest Range and Highest Precision
Measurement Instrument for
Nano-particle Characterization "nano partica SZ-100V2 series"*

Nanoparticle Analyzer
nano partica SZ-100V2 series



A highly advanced analyzer solves the mysteries of the nano-world. A single device analyzes the three parameters that characterize nanoparticles: particle size, zeta potential, and molecular weight.

CE

CE marking compliant

Newly developed to satisfy the need for devices to simply and accurately evaluate the size and dispersion stability of nanoparticles, the key to nanotechnology advancement:

nano partica **SZ-100V2** Series Nanoparticle Analyzer

Nanotechnology research and development is a continuously evolving effort to control substances at the atomic and molecular level in order to achieve new and better materials and products.

The miniaturization of components – that is, control at the nanolevel – is necessary to achieve faster, higher-performance devices and functions and to reduce energy consumption. Nanotechnology has come to play a key role in wide-ranging fields that affect our daily lives, including food, cosmetics, and the life sciences.

Clear and simple multi-parameter analysis of nanoparticles!
Three analyzers in a single compact body deliver high-sensitivity,
high-accuracy analysis of each measurement parameter.



Particle Size Measurement Range $0.3 \text{ nm to } 10 \mu\text{m}$

The SZ-100V2 Series measures particle size and particle distribution width by dynamic light scattering (DLS).

Analysis across a wide range of sample concentrations: Measurement of samples ranging from low ppm-order concentrations to high-concentration samples in double-digit percentages is possible. Accepts commercially available sampling cells. Analysis of small-volume samples is also possible.



Zeta Potential Measurement $-500 \text{ to } +500 \text{ mV}$

Analysis of sample volumes as small as 100 μL using HORIBA-developed microelectrophoresis cells. Use the value of zeta potential to predict and control dispersion stability. High zeta potential magnitudes indicate a stable dispersion, useful for formulation work.



Molecular weight $1 \times 10^3 \text{ to } 2 \times 10^7 \text{ Da}$

Absolute molecular weight (Mw) and the second virial coefficient (A_2) are obtained by performing static light scattering measurement as a function of sample concentration and preparing Debye plots.

The SZ-100V2 Series applies sophisticated intelligence and learning capability to rapidly determine nanoparticle properties!

- Since the SZ-100V2 Series analyzer covers a wide sample concentration measurement range, sample dilution and other preprocessing is nearly eliminated. The use of a dual optical system enables measurement of high-concentration samples such as slurry and ink pigments as well as low-concentration proteins and polymers.
- A single device analyzes the three parameters that characterize nanoparticles: particle size, zeta potential, and molecular weight.
- HORIBA-developed disposable cells for zeta potential measurement prevent sample contamination. Simple analysis by means of ultra micro-volume dedicated cells (volume as low as 100 μL). Suitable for analysis of dilute samples.
- HORIBA-developed electrode for zeta potential cell made from carbon material, the material is not corroded by high salt samples such as saline.



Simple and Convenient Operation

Simply fill the sample cell and place the cell in the analyzer.

A space-saving body design makes the analyzer suitable for installation in any laboratory environment.

Operation Procedure

01



[Sampling]

Fill the sample cell.

02



[Cell Set-Up]

Place the cell in the analyzer.

03



[Start Measurement]

Click the Start button.

04



[Results Display]

The measurement results are displayed.




Maintenance-Free

No maintenance or cleaning of the analyzer is required. After measurement, simply clean or dispose of the cell.




Particle Size Measurement Principle


The SZ-100V2 uses the technique of dynamic light scattering to determine particle size. Dynamic light scattering is the measurement of fluctuations in scattered light intensity with time. These fluctuations in intensity arise due to the random Brownian motion of the nanoparticles. Therefore, the statistical behavior of these fluctuations in scattered intensity can be related to the diffusion of the particles. Since larger particles diffuse more slowly than small particles one can readily relate particle size to measured fluctuations in light scattering intensity. With modern instruments such as the SZ-100V2 the technique is rapid and reliable.

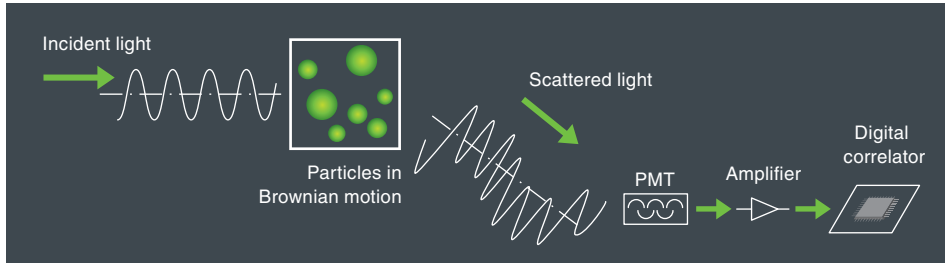

Graphic Rendering of Diffusing Nanoparticles

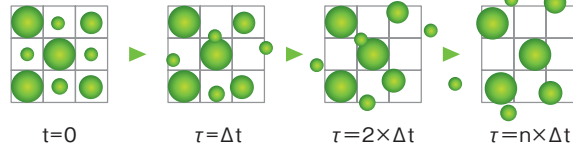

The Relationship Between the Autocorrelation Function and Particle Size

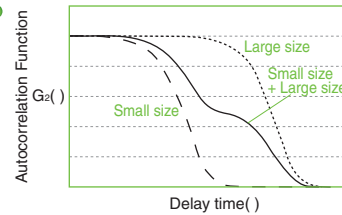

Measurement of the autocorrelation function is done by comparing the scattered light intensity at some reference time t and after some delay time τ . For a very short delay time, the particles have not had a chance to move and therefore the scattered light intensity is unlikely to change much. So, the autocorrelation function has a high value. For a very long delay time, the particles have had a chance to move significantly, and the autocorrelation function has a low value. This low value is related to the time average scattered intensity. The rapidity of this decay from high values to low values corresponds to the speed of particle motion and therefore to the particle size. The measured autocorrelation function typically has an exponential decay and the diffusion coefficient can be calculated with the following (simplified) relationship


$$G^{(2)}(\tau) = B + Bf \exp(-2D_m q^2 \tau)$$

$G^{(2)}(\tau)$: Measured amplitude autocorrelation function B : So-called baseline f : Instrument constant D_m : Particle diffusion coefficient q : Scattering vector given by $(4\pi n/\lambda)\sin(\theta/2)$ τ : Delay time

Particle size is calculated from the diffusion coefficient using the Stokes-Einstein equation:

$$D_h = kT / (3\pi\eta D_m)$$

D_h : Particle hydrodynamic size k : Boltzman constant T : Thermodynamic temperature and η : Viscosity

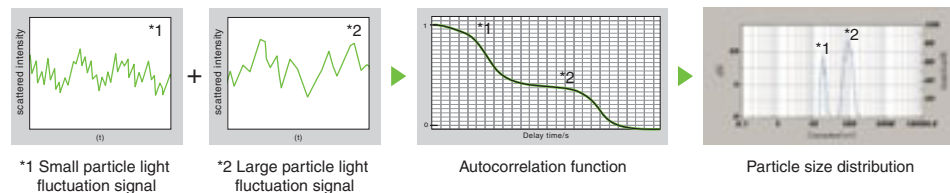

Features of HORIBA's Optical System

1 High Sensitivity Optical Components

The key to accurately and rapidly evaluating size with dynamic light scattering is to use a high-energy laser light source and a sensitive detector. HORIBA uses a green laser. Scattering intensity is inversely proportional to the fourth power of wavelength. Therefore, the green laser gives more scattering intensity per milliwatt than the more commonly used red laser. Since avalanche photodiodes, APD's, are less sensitive to green light and photomultiplier tubes PMT's, are more sensitive to green light, HORIBA has included the most sensitive PMT detector available. In addition, the dead time of a PMT is shorter than that of an APD and therefore the PMT detector dynamic range is superior.


2 Conformance with Standards

The SZ-100V2 series conforms to ISO 13321:1996 and JIS Z8826:2005.

3 Automatic Measurement Optimization

The analyzer features the ability to measure particle size under a number of conditions. In order to eliminate guesswork, measurement conditions can be automatically selected for each sample by using data obtained from that sample.


The three angle system of the SZ-100V2 enables analysis of a wide range of high concentration and dilute samples

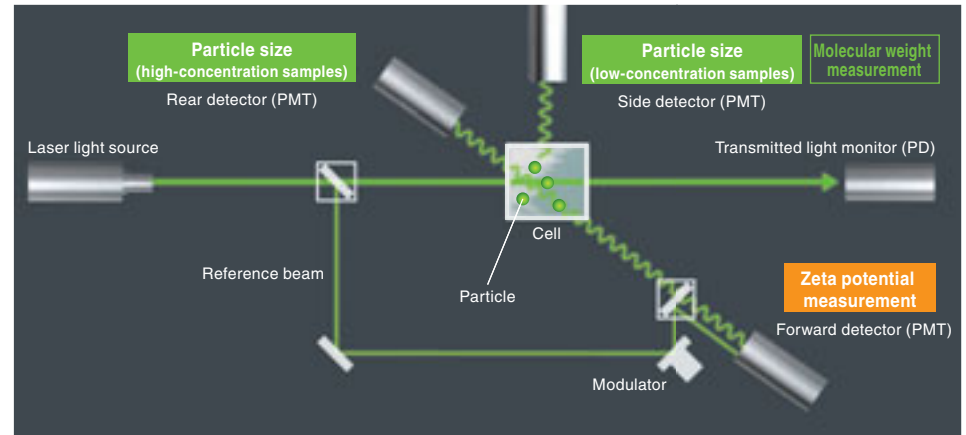

The optical configurations shown can be automatically selected depending on the concentration of the analyzed sample.


High Concentration Samples

In order to minimize the effect of multiple scattering the analyzer detects back-scattered light from a scattering volume close to the cell wall.

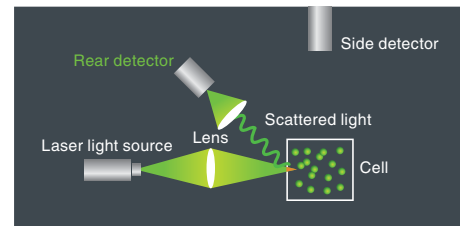

Dilute Samples

In order to minimize the effect of stray light and maximize signal to noise ratio the analyzer detects scattered light at a right angle.

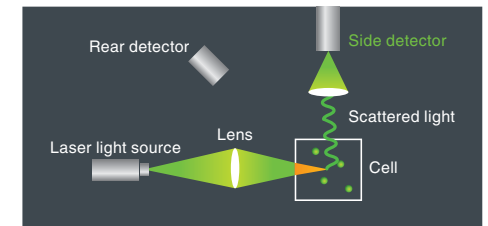




Zeta Potential Measurement Principle

(Laser Doppler Electrophoresis)

Many nanoparticles or colloidal particles have a surface charge when they are in suspension. When an electric field is applied, the particles move due to the interaction between the charged particle and the applied field. The direction and velocity of the motion is a function of particle charge, the suspending medium, and the electric field strength. Particle velocity is then measured by observing the Doppler shift in the scattered light. The particle velocity is proportional to the electrical potential of the particle at the shear plane which is the zeta potential. Thus, this optical measurement of particle motion under an applied field can be used to determine zeta potential.

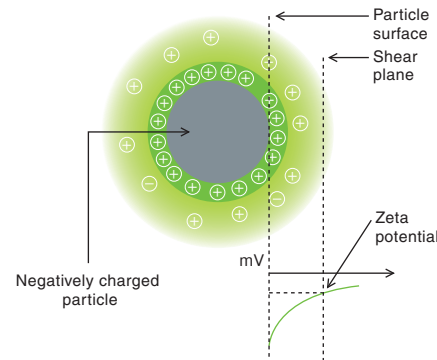

Electrophoresis


Particle motion under an applied electric field is known as electrophoresis. The method used by the SZ-100V2 is known as laser Doppler electrophoresis. Sample particles are suspended in a solvent of known refractive index, n , viscosity, η , and dielectric constant, ϵ . The sample is irradiated with laser light of wavelength λ . An electric field with strength E is applied. Due to the electric field, the particles are moving. Since the particles are moving, the scattered light has a frequency (Doppler) shift proportional to the particle charge. The frequency shift of the scattered light at angle θ is measured and the particle velocity V is determined from the frequency shift. Mobility is then readily obtained as the ratio of velocity to electric field strength V/E . Zeta potential is then found from mobility using a model, the most common of which is the Smoluchowski model.

$$U = \frac{\lambda v a}{2E n \sin(\theta/2)}$$

The following equation is used for the relationship between the calculated electrical mobility and zeta potential.

$$\zeta = \frac{U \eta}{\epsilon f(\kappa a)}$$

ζ : Zeta potential U : Electrical mobility E : Electric field strength n : Solvent refraction index ϵ : Solvent dielectric constant
 η : Solvent viscosity $f(\kappa a)$: Henry coefficient

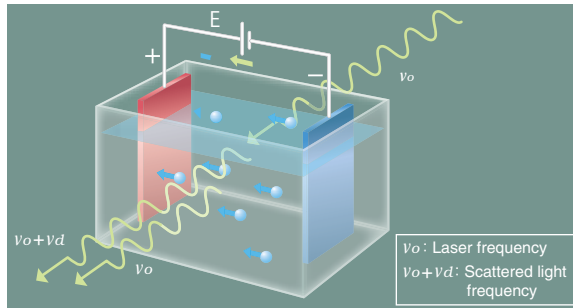


Features

- 1 Extremely low sample volume makes it possible to measure precious or rare samples.
- 2 Modern signal processing electronics efficiently convert optical signals to mobility and zeta potential information. There is no need to manually calculate particle velocity or match speeds.

3 Cell design minimizes electro-osmotic flow to enhance sensitivity.

Particles are not the only objects that acquire a surface charge when in contact with a liquid. Macroscopic objects such as cell or capillary walls do as well. Due to electrostatic attraction, ions with a charge opposite to that of the wall will accumulate close to the wall. And, when an electric field is applied during zeta potential measurement these ions will move in response to the applied field. The moving ions drag the fluid along, creating bulk flow called electro-osmotic flow. This flow will disturb particle motion and distort zeta potential measurements. By eliminating the capillary between the electrodes, the HORIBA zeta potential cell minimizes this effect and maximizes instrument sensitivity.

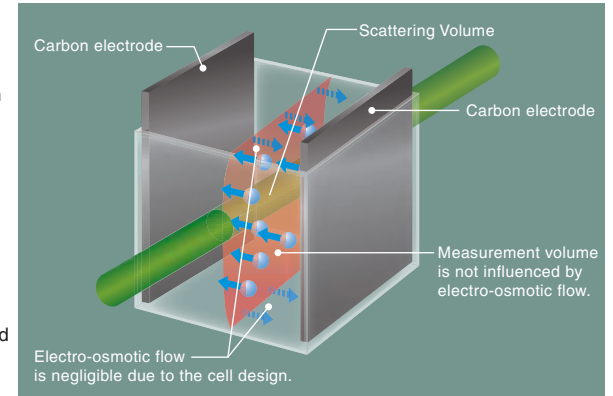


Molecular Weight Measurement Principle

Molecular weight of macromolecules such as polymers, proteins, or starches is determined in two ways with the SZ-100V2. The first method is the use of the dynamic light scattering size information and the empirical Mark Houwink Sakurada equation. The second method is analysis with a Debye plot. Both of these methods are described below.

The Mark Houwink Sakurada equation relates the diffusion coefficient determined by dynamic light scattering to the molecular weight. All that is required are two empirical constants for the selected polymer-solvent system, an exponent and a prefactor. If the constants are not in the SZ-100V2 software database, the user can add new constants for rapid analysis. This technique has the advantage that sample concentration need not be well known.


The Debye plot is obtained by first measuring the excess static light scattering intensity of a series of solutions with well known concentration. Here, the excess intensity refers to the increase of the scattered intensity of the solution compared to the pure solvent. Plotting a quantity proportional to the concentration over the excess scattering as a function of concentration yields a straight line. Extrapolating to zero concentration yields the reciprocal of molecular weight. The graph below shows a typical result.

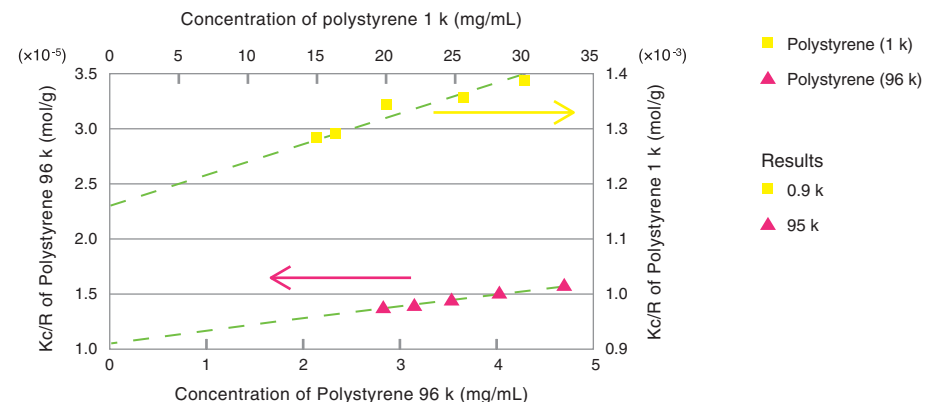

The operator selects a measurement mode (particle size, zeta potential, or molecular weight), loads the sample when the measurement screen appears, and begins measurement. The SZ-100V2 Series offers the ultimate in clear, simple operability. 21 CFR Part 11 software is available.

● Quick and Simple Operation

Measurement conditions are readily set manually or with user programmable methods that can be tied to custom buttons. Operators need merely click a button to begin.

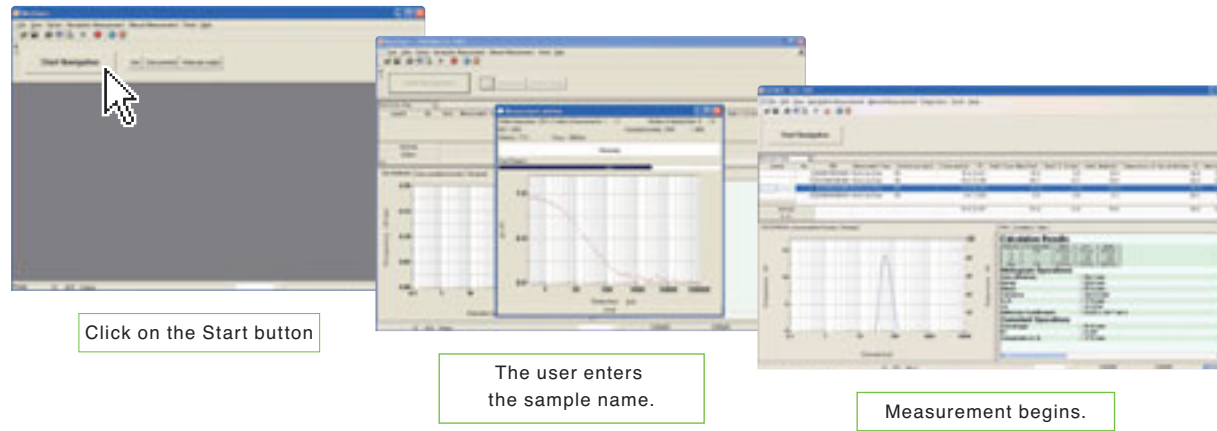


Click on the Start button

The user enters the sample name.

Measurement begins.

● Navigation Creation Is Simple

Use the software wizard to select analysis conditions. If desired, assign a button for fast analysis in the future.

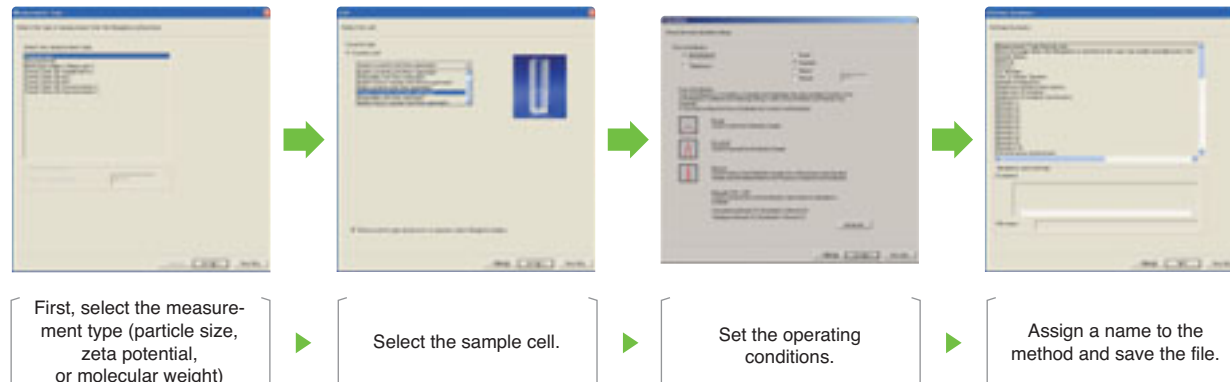

First, select the measurement type (particle size, zeta potential, or molecular weight)

Select the sample cell.

Set the operating conditions.

Assign a name to the method and save the file.

The software follows a progression of selecting measurement conditions and procedures and creates a navigation file.

● Measurement Accuracy

HORIBA confirms measurement performance prior to product shipment using HORIBA-approved standard samples to confirm accuracy and reproducibility as per the tables below. To ensure high-level, stable performance, HORIBA delivers products manufactured in accordance with rigorous quality control systems worldwide.

Particle size

Particle size measurement accuracy using NIST-traceable polystyrene latex standards particles is as shown below.

Particle size standard value (nm)	Concentration	Standard
100 nm	100 ppm	Measured values for cumulant average size are within $\pm 2\%$. (This does not include variation in the standard particles themselves.)

Particle size measurement reproducibility is as shown below.

Particle size standard value (nm)	Concentration	Standard
100 nm	100 ppm	The CV value for 6 repeated measurements is less than 2%.
100 nm	10 wt %	The CV value for 6 repeated measurements is 5% or less.

*Conforms to ISO 13321: 1996, ISO 22412 : 2008 and JIS Z 8826: 2005.

Zeta Potential

Using a HORIBA-designated colloidal silica sample, HORIBA confirms that the measured value is higher than -75 mV and lower than -40 mV. Reproducibility for 6 repeated measurements is within 10% or less in CV value.


Molecular Weight

The measured value is within $\pm 10\%$ of the standard value using a polystyrene standard sample (Nominal molecular weight: 96,000).


Applications

● Biomaterials: Gold colloid particle size measurement results


Au colloids (NIST)	RM8011	RM8012	RM8013
Nominal Size (nm)	10	30	60
NIST reference size by dynamic light scattering (nm)	13.5	26.5	55.3
Size measured with SZ-100V2 (nm)	11.0	26.6	55.4

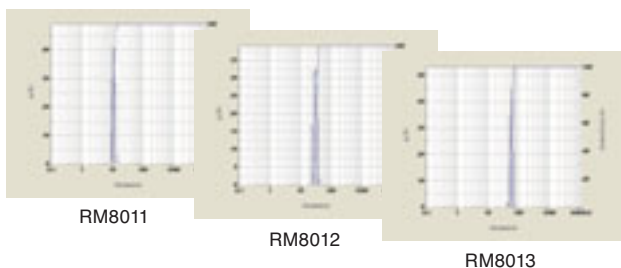

● Gold colloid particle (2 nm) size measurement results (with high power laser 532 nm 100 mW)

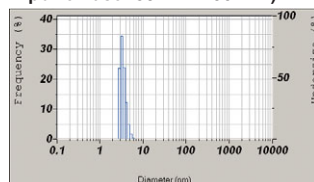


● Lysozyme (from egg white) particle size measurement result (with high power laser 532 nm 100 mW)

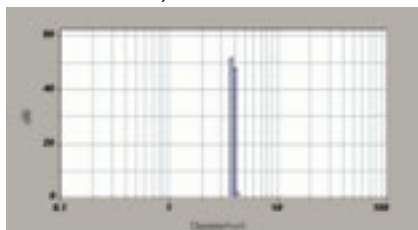


Sample concentration: 0.05 mg/mL
Acetic acid buffer: pH = 4.3
Average diameter: 4.0 nm

● Thiamin hydrochloride (Vitamin B1 hydrochloride) particle diameter measurement result

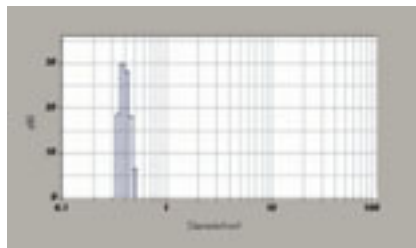


Sample concentration: 300 mg/mL
Average diameter: 0.4 nm

● NIST SRM 1980 -FeOOH zeta potential measurement result

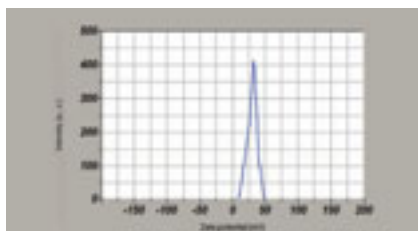

Sample concentration: 50 ppm, pH = 2.5
Mobility (rated): $2.53 \pm 0.12 \mu\text{m} \cdot \text{cm}/\text{Vs}$
Measurement results: Mobility = $2.53 \mu\text{m} \cdot \text{cm}/\text{Vs}$
Zeta potential = 32.9 mV

● Isoelectric point of silica measurement result

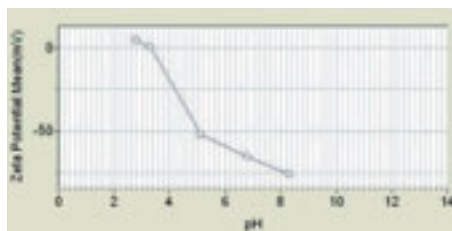


Sample concentration: 0.01 mol/L
(Adjusted to 10 w% with KCl)
Zeta potential = -38.3 mV

Accessories

● Sample Cell Types and Specifications

We can guide you in selecting the right cell for your application.

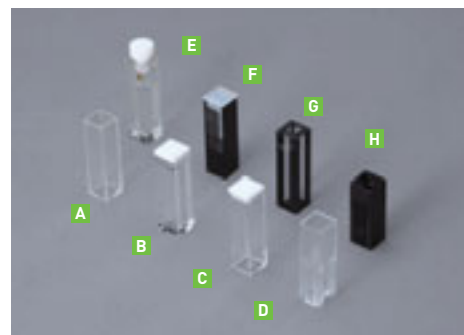

Zeta potential measurement disposable cells
(For zeta potential and particle size measurement, 100 μL , Aqueous)

	Cell Name	Measurement Application	Remarks
A	Disposable cell	Particle size/ Molecular weight	Plastic, 4 surfaces clear, 100 pieces, Full volume 4000 μL (Minimum sample volume 1000 μL)
B	Semi-micro cell	Particle size	Quartz, 4 surfaces clear, Full volume 1600 μL (Minimum sample volume 400 μL)
C	Glass cell	Particle size/ Molecular weight	Glass, 4 surfaces clear, Full volume 4000 μL (Minimum sample volume 1000 μL)
D	Semi-micro disposable cell	Particle size	Plastic, 2 surfaces clear, 100 pieces, Full volume 800 μL (Minimum sample volume 400 μL)
E	Cell with lid	Particle size/ Molecular weight	Quartz, 4 surfaces clear, Full volume 4000 μL (Minimum sample volume 1000 μL)
F	Micro-cell (Side detector only)	Particle size/ Molecular weight	Quartz, 3 surfaces clear, Side detector only, Full volume 30 μL (Minimum sampling volume 12 μL)
G	Sub-micro cell	Particle size/ Molecular weight	Quartz, 4 surfaces clear, Full volume 750 μL (Minimum sampling volume 250 μL)
H	Flow cell	Particle size/ Molecular weight	Quartz, 3 surfaces clear, Full volume 100 μL (Minimum sampling volume 100 μL), 2 connectors with pH controller
I	Zeta potential plastic cell	Zeta potential	For aqueous sample, 20 pieces
J	Zeta potential glass cell	Zeta potential	For organic solvent, 50 replacement gold electrodes, PTFE lid, and 2 caps.

● Autotitrator

This device can be used to automatically prepare plots of zeta potential or particle size as a function of pH. It is an excellent choice for iso-electric point determination.



pH Controller Accessory Specifications

- Number of titrant bottles: 2
 - pH adjustment range: 1 - 13
 - Burette volume: 5 mL
 - Power supply: AC 100 - 240 V, 50/60 Hz, 45 VA
 - Dimensions and weight:
Body: 141 (W) x 296 (D) x 367 (H) mm, approx. 4 kg
Circulation Pump: 124 (W) x 202 (D) x 122 (H) mm, approx. 1.7 kg
 - Option: Magnetic Stirrer
- pH electrode part number 3200585428
pH calibration unit parts number 3200043642

SZ-100-S2 Measurement Specifications

Model	SZ-100-S2 and SZ-100-HS2 (particle size and molecular weight measurement only)
Measurement principles	Particle size measurement: Dynamic Light Scattering Molecular weight measurement: Debye plot method (static scattered light intensity)
Measurement range	Particle size: 0.3 nm to 10 μm Molecular weight: 1000 to 2 x 10 ⁷ Da (Debye plot) 540 to 2 x 10 ⁷ Da (MHS Equation) ^{*1}
Maximum sample concentration	40 wt% ^{*2}
Particle size measurement accuracy	Measurement accuracy of ±2% for NIST traceable polystyrene latex 100 nm spheres (not including variation in the standard particles themselves)
Measurement angles	90° and 173° (automatic or manual selection)
Cells	Cuvettes
Measurement time	Approx. 2 min. under ordinary conditions (from the start of measurement to the display of results for particle size measurement)
Required sample volume	Minimum volume of 12 μL ^{*3} to 1000 μL (differs depending on cell material)
Usable liquids	Water, ethanol, organic solvents

*1: Mark-Howink-Sakurada Equation, depending on sample. *2: Depending on sample. *3: F Micro-cell.

SZ-100-Z2 Measurement Specifications

(Particle size and molecular weight measurement specifications are the same as for the SZ-100-S2.)

Model	SZ-100-Z2 and SZ-100-HZ2 (with zeta potential measurement unit)
Measurement principle	Zeta potential measurement: Laser Doppler electrophoresis
Measurement range	-500 to +500 mV
Size range suitable for measurement	Minimum 2.0 nm, Maximum 100 μm ^{*4}
Measurement conductivity range	0 to 20 S/m ^{*5}
Maximum sample concentration	40 wt% ^{*6}
Cells	Dedicated cell with electrodes
Measurement time	Approx. 2 min. under ordinary conditions
Required sample volume	100 μL
Carrier fluids	Water

*4: Depending on sample. *5: Recommended sample conductivity range: 0 to 2 S/m. *6: Depending on sample.

Analyzer Specifications (SZ-100-S2, SZ-100-Z2, SZ-100-HS2 and SZ-100-HZ2)


Measuring unit optical system	Light source: Diode pumped frequency doubled laser (532 nm, S2 / Z2 10 mW, HS2 / HZ2 100 mW) Detectors: Photomultiplier tubes (PMT)
Laser classification	Class I
Operating temperature and humidity	15 - 35 °C, RH 85 % or less (no noncondensing)
Holder temperature control temperature settings	0 - 90 °C (up to 70 °C for cells with electrodes and plastic cells)
Purging	Dry gas purge port tube connection is possible.
Power supply	AC 100 - 240 V, 50/60 Hz, 150 VA
Dimensions	528 (W) x 385 (D) x 273 (H) mm (excluding protrusions)
Weight	25 kg
Personal computer	Windows computer with one available USB port
Interface	USB 2.0 (between measuring unit and PC)
OS	Windows® 10 32/64 bit

Class I laser product CE certification



* Composite photographs are inserted into the PC screens.

Dimensions (mm)

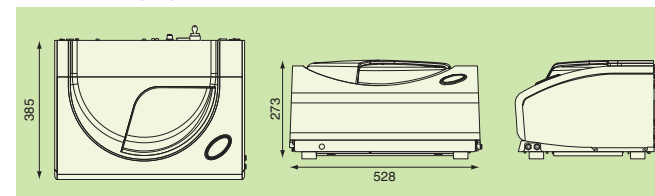

Data Processing

Navigation files turn complex parameter input into simple to use operating procedures. / Store 100 data items on a data list. / Display individual data items with a single mouse click. / Perform pH, temperature, and sample concentration trend measurement.

Particle Size Measurement

Real-time display of the autocorrelation function / Display of median size, specific surface area, mode size, average size, standard deviation, coefficient of variation, span value, percentage size (max. of 10 items), Z average, polydispersity index, size percentage (max. of 10 items displayed) / Particle distribution graph, autocorrelation function, residual error / Refractive index, viscosity, computing range, and data recalculation after measurement

Molecular Weight Measurement

Molecular Weight Measurement
Real-time display of Debye plots / Display of molecular weight and the second virial coefficient / Recalculation of Debye plot graph display data

Zeta Potential Measurement

Zeta potential, standard deviation, electrophoretic mobility, and average zeta potential at each peak / Display of zeta potential graphs, mobility graphs, recalculation of data

Options

21CFR Part 11 software / Zeta potential measurement organic solvent cells / pH control unit / IQ/OQ/PQ support



The HORIBA Group adopts IMS (Integrated Management System) which integrates Quality Management System ISO9001, Environmental Management System ISO14001, and Occupational Health and Safety Management System OHSAS18001. We have now integrated Business Continuity Management System ISO22301 in order to provide our products and services in a stable manner, even in emergencies.



Please read the operation manual before using this product to assure safe and proper handling of the product.

- The specifications, appearance or other aspects of products in this catalog are subject to change without notice. ●Please contact us with enquiries concerning further details on the products in this catalog.
- The color of the actual products may differ from the color pictured in this catalog due to printing limitations. ●It is strictly forbidden to copy the content of this catalog in part or in full. ●The screen displays shown on products in this catalog have been inserted into the photographs through compositing.
- All brand names, product names and service names in this catalog are trademarks or registered trademarks of their respective companies.

<http://www.horiba.com>

HORIBA, Ltd. Japan

Head Office
2 Miyahogashi, Kisshoin, Minami-ku, Kyoto, Japan
Phone: 81 (75) 313-8123 Fax: 81 (75) 321-5725

HORIBA KOREA Ltd. Korea

25, 94-Gil, Ijjik-Ro, Manan-Gu, Anyang-Si, Gyeonggi-Do, 13901, Korea
Phone: 82 (31) 296-7911 Fax: 82 (31) 296-7913

HORIBA India Private Limited India

246, Okhla Industrial Estate, Phase 3 New Delhi-110020, India
Phone: 91 (11) 4646-5000 Fax: 91 (11) 4646-5020
Bangalore Office
No.55, 12th Main, Behind BDA Complex, 6th sector, HSR Layout, Bangalore South, Bangalore-560102, India
Phone: 91 (80) 4127-3637

HORIBA Taiwan, Inc. Taiwan

8F.-8, No.38, Taiyuan St. Zhubei City, Hsinchu County 30265, Taiwan (R.O.C.)
Phone: 886 (3) 560-0606 Fax: 886 (3) 560-0550

HORIBA (China) Trading Co., Ltd. China

Unit D, 1F, Building A, Synnex International Park, 1068 West Tianshan Road, Shanghai, 200335, China
Phone: 86 (21) 6289-6060 Fax: 86 (21) 6289-5553
Beijing Branch
12F, Metropolis Tower, No.2, Haidian Dong 3 Street, Beijing, 100080, China
Phone: 86 (10) 8567-9966 Fax: 86 (10) 8567-9066
Guangzhou Branch
Room 1611 / 1612, Goldlion Digital Network Center, 138 Tiyu Road East, Guangzhou, 510620, China
Phone: 86 (20) 3878-1883 Fax: 86 (20) 3878-1810

HORIBA (Thailand) Limited Thailand

393, 395, 397, 399, 401, 403 dad Ya Road, SomdetChaophraya. Khlongsan District Bangkok. Khlongsan District 10600 Thailand.
Phone: 66 (2) 861-59-95 ext.123 Fax: 66 (2) 861-52-00
East Office
850 / 7 Soi Lat Krabang 30 / 5, Lat Krabang Road, Lat Krabang, Bangkok 10520, Thailand
Phone: 66 (2) 734-4434 Fax: 66 (2) 734-4438

HORIBA Instruments Incorporated USA

9755 Research Drive, Irvine, CA 92618, U.S.A.
Phone: 1 (949) 250-4811 Fax: 1 (949) 250-0924
Piscataway Office
20 Knightsbridge Road, Piscataway, NJ 08854, U.S.A.
Phone: 1 (732) 494-8660 Fax: 1 (732) 549-5125

HORIBA Instruments Brasil, Ltda. Brazil

Rua Presbitero Plinio Alves de Souza, 645, Parte A, Loteamento Multivias, Jardim Ermida II - Jundiá Sao Paulo - CEP 13.212-181 Brazil
Phone: 55 (11) 2923-5400 Fax: 55 (11) 2923-5490

HORIBA Vietnam Co., Ltd. Vietnam

Lot 3 and 4, 16 Floor, Detech Tower II, No. 107 Nguyen Phong Sac Street, Dich Vong Hau Ward, Cau Giay District, Hanoi, Vietnam
Phone: 84 (24) 3795-8525 Fax: 84 (24) 3795-8553

PT HORIBA Indonesia Indonesia

Jl. Jalur Sutera Blok 20A, No.16-17, Kel. Kunciran, Kec. Pinang Tangerang-15144, Indonesia
Phone: 62 (21) 3044-8525 Fax: 62 (21) 3044-8521

HORIBA FRANCE SAS France

16-18, rue du Canal, 91165, Longjumeau Cedex, France
Phone: 33 (1) 69-74-72-00 Fax: 33 (1) 69-09-07-21

HORIBA Jobin Yvon GmbH Germany

Neuhofstrasse 9, D_64625, Bensheim, Germany
Phone: 49 (89) 62-51-84-750 Fax: 49 (89) 62-51-84-7520

HORIBA ITALIA SRL Italy

Via Luca Gaurico 209-00143, Roma, Italy
Phone: 39 (6) 51-59-22-1 Fax: 39 (6) 51-96-43-34

HORIBA UK Limited UK

Northampton Office
Kyoto Close Moulton Park, Northampton NN3 6FL, UK
Phone: 44 (1604) 542-500 Fax: 44 (1604) 542-699

HORIBA Instruments (Singapore) Pte Ltd. Singapore

3 Changi Business Park Vista #01-01, Akzonobel House, Singapore 486051
Phone: 65 (6) 745-8300 Fax: 65 (6) 745-8155

Bulletin:HRE-3677F

Printed in Japan 1912SK152