

Application News

Ion-Pair Reversed-Phase LC/MS Analysis of siRNA under Denaturing and Non-Denaturing Conditions

Junna Nakazono

User Benefits

- ◆ The LCMS-9050 quadrupole time-of-flight mass spectrometer can be used to characterize siRNA.
- ◆ The elution of siRNA in double-stranded (non-denatured) or single-stranded (denatured) state can be confirmed by changing the column temperature.
- ◆ The LabSolutions Insight™ Biologics analysis software enables analysis of multiple oligonucleotide sequences at once.

Introduction

Among oligonucleotide therapeutics, antisense and siRNA therapeutics are actively studied as new modalities for treating genetic and intractable diseases. Unlike single-stranded antisense, siRNA is double-stranded RNA and known to separate into single-stranded RNAs under heated or alkaline conditions. Quality control of siRNA drugs requires characterization of both single-stranded and double-stranded RNA.

This application describes an example of LC/MS analysis of siRNA under denaturing and non-denaturing conditions in ion-pair reversed-phase chromatography. Multiple nucleic acid sequences (sense and antisense strands) were analyzed at the same time using LabSolutions Insight Biologics, a dedicated software for oligonucleotide characterization.

Samples

Double-stranded siRNA annealed with sense and antisense strands of the following sequences was used.

Sense:

G-Um-A-A-Cm-Cm-A-A-G-A-G-Um-A-Um-Um-Cm-Cm-A-Um-dT-dT

Antisense:

A-U-G-G-A-A-Um-A-C-U-C-U-U-G-G-U-Um-A-C-dT-dT

(d: 2'-deoxy and m: 5-methyl)

Analytical Conditions

Analysis was performed with Nexera™ XS inert UHPLC and LCMS-9050 quadrupole time-of-flight mass spectrometer systems. The analytical conditions are shown in Table 1.

Table 1 Analytical Conditions

UHPLC (Nexera XS inert)	
Column:	Shim-pack Scepter Claris C18-300*1 (100 mm x 2.1 mm I.D., 1.9 μm)
Mobile Phase A:	100 mM HFIP, 10 mM TEA - water
Mobile Phase B:	100 mM HFIP, 10 mM TEA - methanol
Gradient Program:	B Conc. 5 % (0 min) – 50 % (10 min) – 90 % (10.01-12 min) – 5 % (12.1-25 min)
Flowrate:	0.3 mL/min
Column Temp.:	25 or 60 °C
Injection Volume:	1 μL

*1 P/N: 227-31209-02

MS (LCMS-9050)	
Ionization:	ESI negative
Mode:	MS m/z 550-2500
Nebulizing Gas Flow:	2.0 L/min
Drying Gas Flow:	10.0 L/min
Heating Gas Flow:	10.0 L/min
Interface Temp.:	350 °C
DL Temp.:	250 °C
Block Heater Temp.:	400 °C

Configuring the Data Analysis Parameters

LabSolutions Insight Biologics is data analysis software for characterizing oligonucleotides and oligonucleotide impurities. First, the user creates an oligonucleotide sequence in the parameter configuration window using the nucleobases, linkers, ribose and modifications provided by the software. Nucleobases, linkers, ribose, and base modifications can be added and removed in each tab as required. Once an oligonucleotide sequence is entered, the software displays the molecular formula, monoisotopic mass (to the left side), and structural formula (to the right side) of that oligonucleotide (Fig. 1).

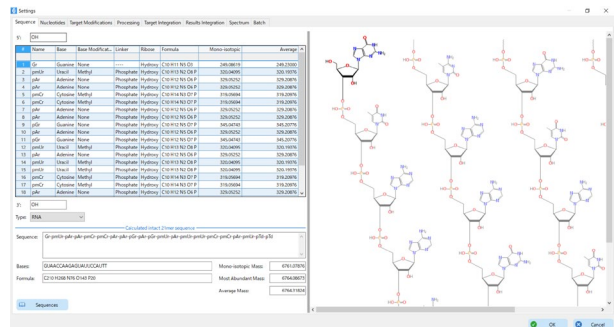


Fig. 1 Parameter Configuration Window

Insight Biologics can analyze multiple sequences. As an example, in Fig. 2, the sequences of sense and antisense strands were set as analysis targets. The specified sequences and optionally added information about nucleobases, linkers, riboses, and base modifications can be saved as an analysis settings file.

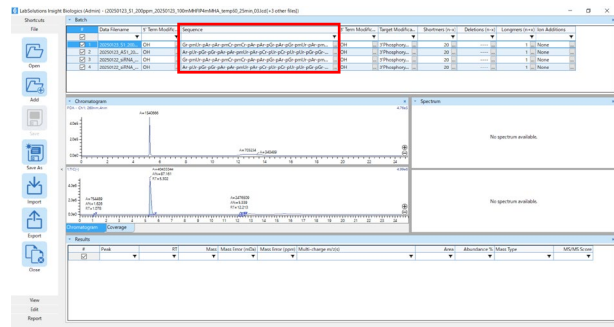


Fig. 2 Example of Multiple Sequence Settings

■ Separation of siRNA under Denaturing and Non-Denaturing Conditions

The UV chromatograms of sense, antisense, and double-stranded siRNA are shown in Fig. 3. The sense and antisense strands were eluted at different retention times at both 60 and 25 °C column oven temperatures. At a 60 °C column oven temperature, double-stranded siRNA was detected as denatured and dissociated into sense and antisense strands. On the other hand, at a 25 °C column temperature, double-stranded siRNA was eluted at a retention time of 7.5 minutes, which was slower than that of sense and antisense strands. In addition, the peak of the sense strand, which seemed to remain without forming a duplex during annealing, was also observed.

■ Results of Identification by LC/MS

Fig. 4 shows the component chromatograms obtained by LC/MS analysis of double-stranded siRNA. In Insight Biologics, the identified oligonucleotide sequence is displayed as a component chromatogram based on the MS1 spectrum and summed with different valences and isotopes. The sense strand and antisense strand sequences were identified at both column temperatures of 60 °C and 25 °C.

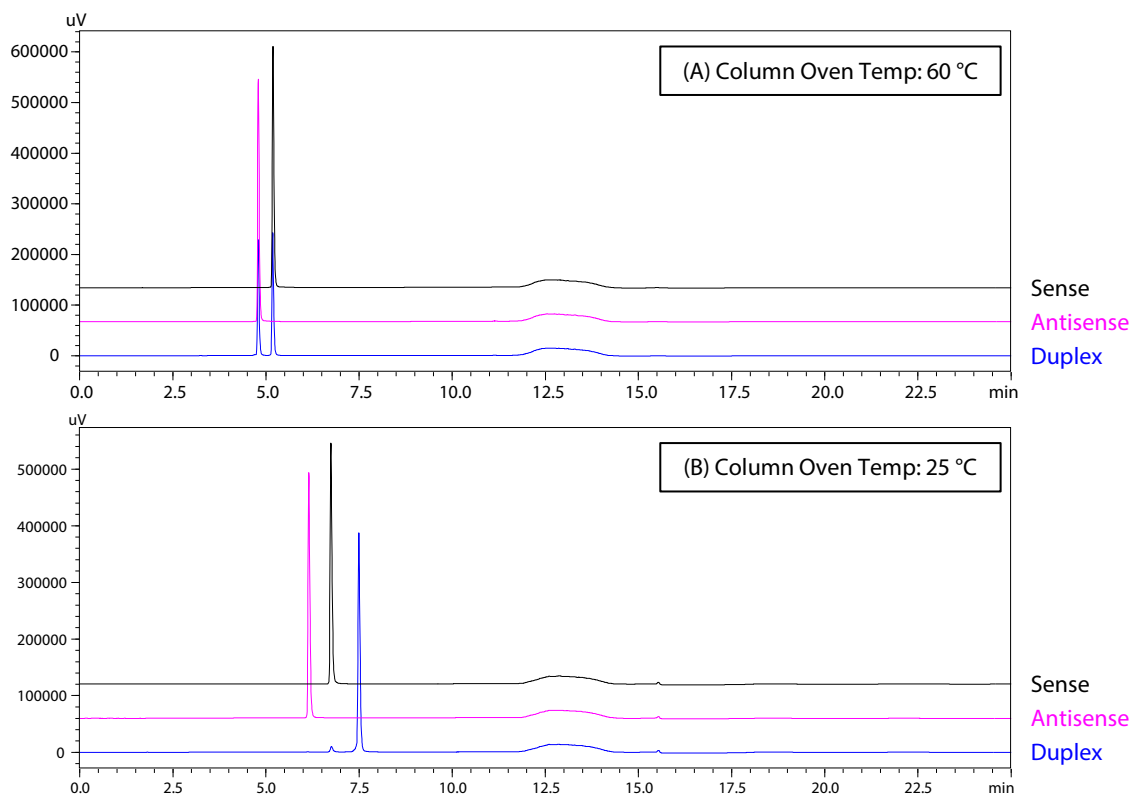


Fig. 3 UV Chromatograms (260 nm) of the Sense Strand, Antisense Strand and Duplex siRNA (A) With 60 °C Column Oven Temp.; (B) With 25 °C Column Oven Temp.

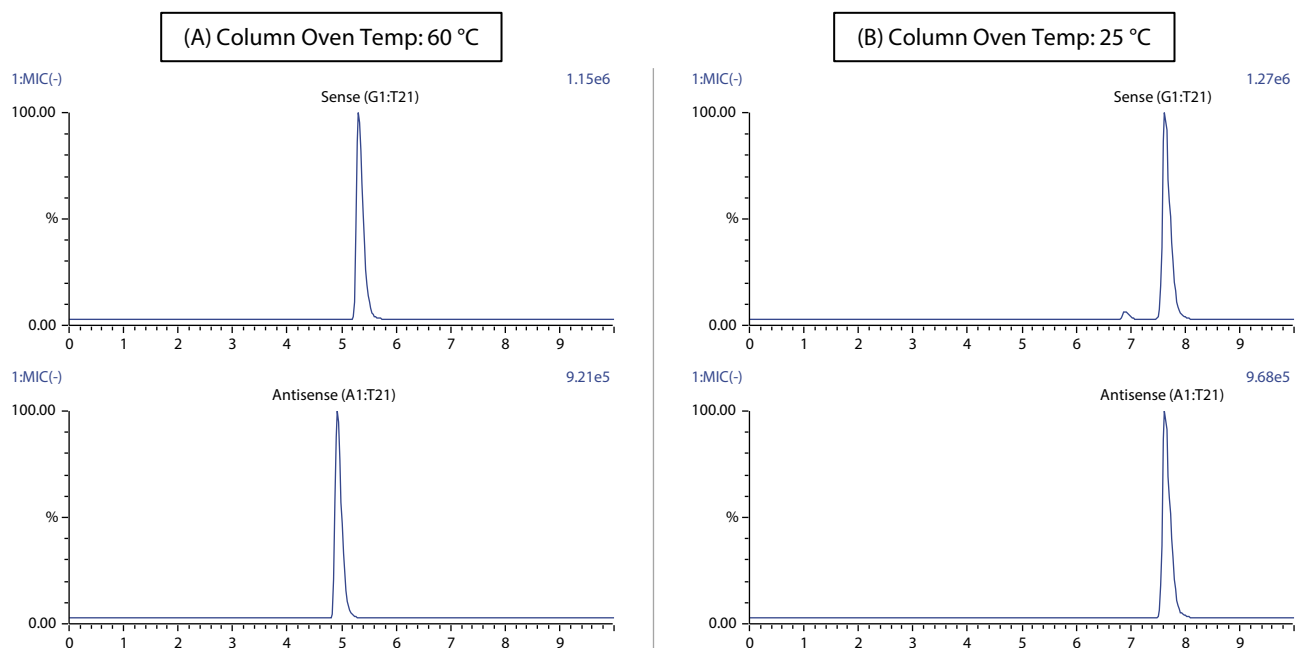


Fig. 4 Component Chromatograms of Duplex siRNA (Results of Identification by LC/MS) (A) With 60 °C Column Oven Temp.; (B) With 25 °C Column Oven Temp.

■ Conclusion

LC/MS analysis of siRNA under denaturing and non-denaturing elution conditions was performed using an LCMS-9050 quadrupole time-of-flight mass spectrometer. Under denaturing conditions at a column temperature of 60 °C, siRNA was eluted in a dissociated single-stranded state, but at a column temperature of 25 °C, siRNA was eluted in a double-stranded state. By using the oligonucleotide analysis software LabSolutions Insight Biologics, multiple nucleic acid sequences can be analyzed at the same time.

■ Acknowledgments

This research was supported by AMED under Grant Number JP21ae0121022, JP21ae0121023, JP21ae0121024 (Project leader: Satoshi Obika).

Related Applications

1. An Oligonucleotide Impurity Analysis Workflow Using LabSolutions Insight™ Biologics Software
[Application News No. 01-00595A-EN](#)
2. Simple Analysis of Impurities in Oligonucleotide Therapeutics Using a Single Quadrupole Mass Spectrometer
[Application News No. 01-00656A-EN](#)
3. Efficient Method Development for Separation of Capped mRNA Fragments
[Application News No. 01-00898-EN](#)

LabSolutions Insight, Nexera, and Shim-pack Scepter are trademarks of Shimadzu Corporation or its affiliated companies in Japan and/or other countries.



SHIMADZU

Shimadzu Corporation

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. See <https://www.shimadzu.com/about/trademarks/index.html> for details.

Third party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

01-00915A-EN

First Edition: Jul. 2025

Revision A: Apr. 2026

[> Please fill out the survey](#)

Related Products

Some products may be updated to newer models.



[> LCMS-9050](#)

Quadrupole Time-of-Flight Liquid Chromatograph Mas...



[> Nexera XS inert](#)

Ultra High Performance Liquid Chromatograph



[> LabSolutions Insight Biologics](#)

Software for Oligonucleotide Sequence Characteriza...

Related Solutions

[> Pharmaceutical and Biopharmaceutical](#)

[> Oligonucleotide and mRNA Therapeutics](#)

[> Price Inquiry](#)

[> Product Inquiry](#)

[> Technical Service / Support Inquiry](#)

[> Other Inquiry](#)